

# SERRES

Helping to cope with the most demanding circumstances, the SERRES Suction System includes special products such as a Specimen Set, Measuring Cup and Splash-Vac. All accessories are compatible with SERRES Suction Bags and are easily connected, disposable and made from environmentally friendly materials.

SERRES Specimen Set is designed for tissue collection. The tissue is collected into a removable sifter inside the cup body. The sifter collects particles larger than 1 mm in size. Samples remain isolated until the laboratory.

SERRES Measuring Cup is designed for the accurate measuring of liquids (e.g. paediatric treatment). Its capacity is 400 ml and measuring tolerance is 5 ml up to 50 ml and 10 ml above. It can be emptied during the procedure without disconnecting the tubes.

SERRES Splash-Vac is ideal for removing large amount of liquids from the floor conveniently. The Splash-Vac can easily be moved either by foot or with a long handle. Its teardrop shape ensures easy access to corners. Splash-Vac is supplied with 3 metres of suction tubing.

## Suction System Accessories




## PRODUCT SPECIFICATION

### Raw materials


Specimen Set:	Polyethylene and polypropylene
Measuring Cup:	SAN, polyethylene and polypropylene
Splash-Vac + tubing:	Polystyrene and polyethylene




## Product information

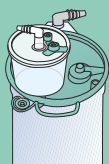

product number	product name	volume	packing size	other information patient port / angle connector
57242	SERRES Specimen Set	300 ml	18 / 108 pcs	12 mm
57247	SERRES Specimen Set	300 ml	18 / 108 pcs	12 / 7,2 mm

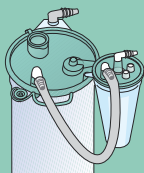

557257	SERRES Measuring Cup	400 ml	12 / 72 pcs	12 / 7,2 mm
--------	----------------------	--------	-------------	-------------

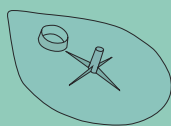

57270	SERRES Splash-Vac + 3 m suction tubing	25 pcs		6 mm
-------	--	--------	--	------